

NEOGENESIS



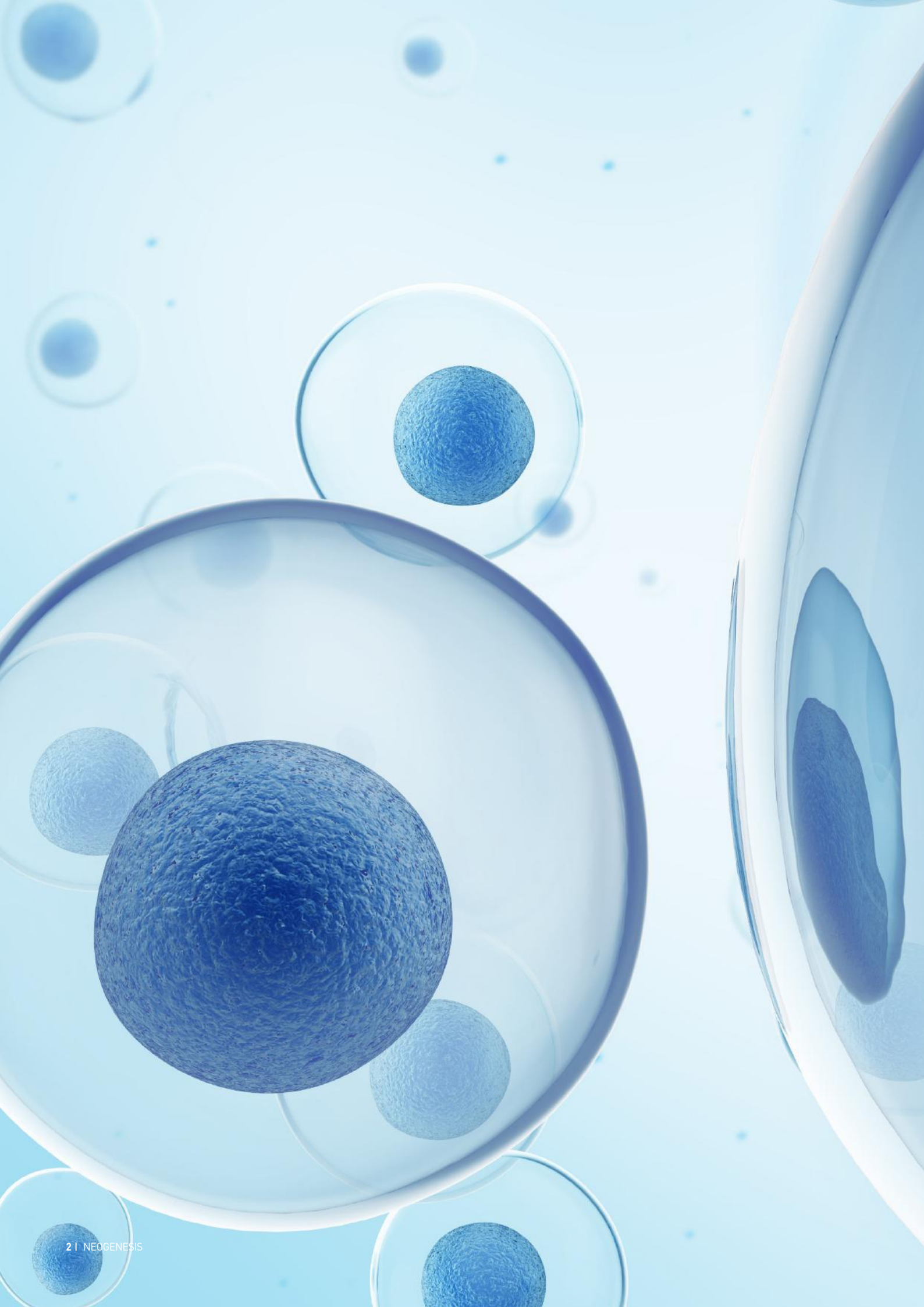


Table of Contents

Company Introduction	4
UNISTATION™ 2ND Ver. - Adipose Derived Stem Cells (ADSCs) Device	6
NeoGenesis PRP Series	10
· Neo PRP (10cc), Genesis PRP (15cc)	
· Pure PRP (20cc), Cure PRP (30cc), Omni PRP (60cc)	
PRP Activator	14
Genesis Centrifuge	15
Plasma Gel - Autologous Filler Making Device	16
REGENOVUE Series - HA Filler	17
· REGENOVUE Series (FINE / FINE+, DEEP / DEEP+, SUB-Q / SUB-Q+)	
REGENOVUE Aqua Shine Series - HA Mesotherapy	19
· REGENOVUE Aqua Shine (Shine / Shine +)	
Neo Filler Cannula	20
Neo Thread Series - PDO Lifting Thread Series	21
· Mono & Screw, Multi, Combi, Cog	
· Nose / Eye Cannula, Double Needle Threads, Molding Thread Cannula	
Beta Series - Skin Care	28
· Beta Scaffold Mask and Beta Serum	
Spotless - Radio Frequency Mole Cleanser Device	30
RIOX - Carboxyl CO₂ Facial Mask	31
Neo - Cain Cream - Lidocaine Cream 10.56%	32

About Us

Company	NeoGenesis Co., Ltd.
CEO	Yongpil, Kim
Establishment	2012
Core Business	R&D, Sales, Export Business of Autologous Cell Therapeutic Products and Medical Devices
Address	Office# 302, 303 E&C Dream Tower, 146, Seonyu-ro, Yeongdeungpo-gu, Seoul, South Korea (07255)
E-mail	info@neogenesis.co.kr
Homepage	www.neogenesis.co.kr
Tel	+82-2-583-1348
Fax	+82-2-2038-0848

NeoGenesis Co., Ltd

We, NeoGenesis Co., Ltd., are known for our professionalism in the autologous cell therapy and stem cell field. We not only strongly focus on our domestic markets, but globally export our medical devices to the Northern/Southern America, Eastern/Central Asia, Europe, and more.

Our mission is to contribute to the healthcare industry for developing autologous cell therapy devices. We continuously research, develop, and extend our knowledge by participating in diverse seminars, exhibitions, trainings, and more.

As we are continuously researching and developing diverse multi-autologous therapeutic stem cell medical devices, we provide safe and effective high-end medical devices for our doctors and patients' high satisfaction and result.

Lastly, we promise to do our best to provide our clients the most safe and professional high-end medical devices and products for a mutual beneficial business relationship.





International Partner Establishment History

2012	2013	2014	2015	2016
United State of America Taiwan	Turkey Netherlands Switzerland Colombia Venezuela Malaysia Thailand China Cyprus United Arab Emirate Kingdom of Saudi Arabia Kuwait Qatar Bahrain Oman Indonesia Czech Republic Chile Jordan Iraq Philippines Iran Russia Turkey	Spain Portugal Greece Italy Russia Estonia Dominica Republic Lithuania Latvia Hong Kong Poland Bulgaria Ukraine Iran	Canada Mexico Iran Vietnam Poland Kuwait Peru Bolivia Ecuador Uzbekistan Australia Crimea India Germany Singapore El Salvador Tajikistan Romania	Moldova Kazakhstan Egypt Azerbaijan Uruguay Slovenia Hungary France Syria Mongolia Bangladesh Denmark



Distributing Over 65 Countries around the Globe






UNISTATION™ 2nd Version

The simplest & easiest way to various cellular therapies

A Unique All-in-One Device

UNISTATION™ is a specially designed medical device to satisfy diverse demands for various autologous cell therapies, including

- Adipose-Derived Stem Cells (ADSCs)
- Stromal Vascular Fraction (SVF)
- Platelet Rich Plasma (PRP)
- Platelet Rich Fibrin (PRF)
- Fat Purification for Lipotransfer

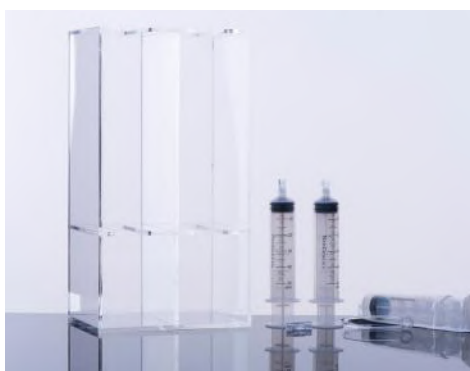

What is Stem Cell?


Stem cells are multi-potent cells and its multi-diversification capability is considered as the core future of medical fields, demanding various cell therapies for a number of incurable diseases.

What is Adipose-Derived Stem Cells (ADSCs)?

Adipose tissue is one of the richest sources of stem cells. ADSCs are actively being researched for various disease treatments.

With assistance of the UNISTATION, ADSCs are easily and abundantly extracted. Furthermore, ADSCs can be used in wide variety of areas as the picture illustrated below.

Application Areas

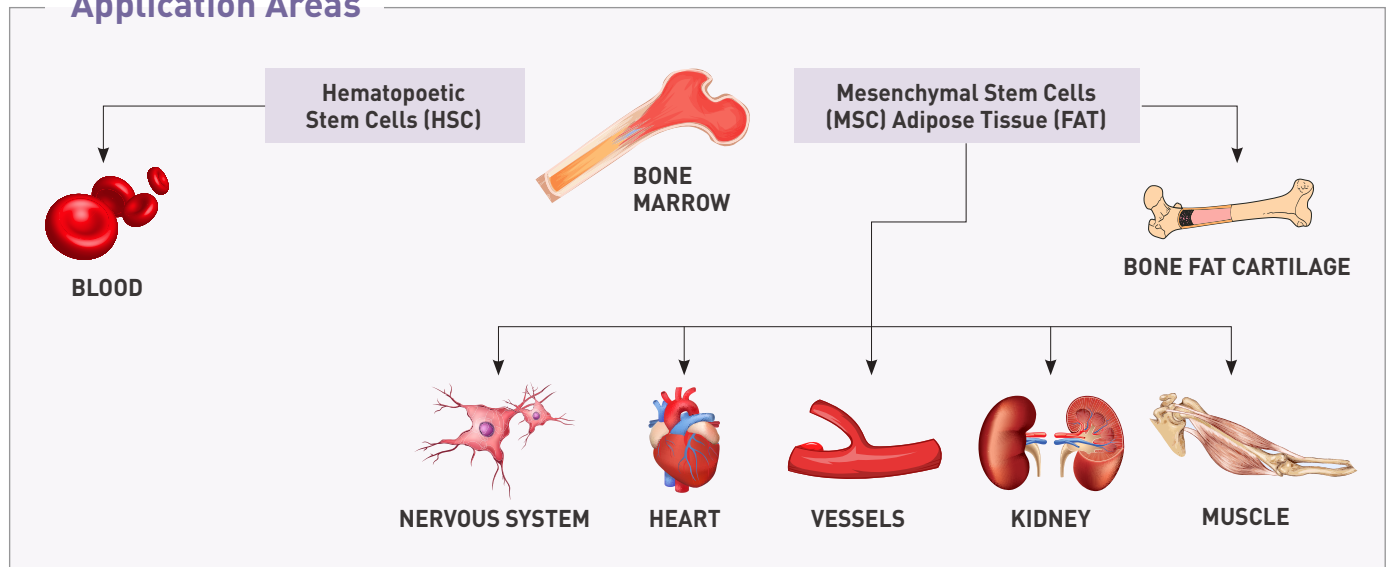

What is Cell-Assisted Lipotransfer (CAL Method)?

CAL Method transfers excessive fat to fill other areas by selecting strong cells among substantial fat cells and processing them through Cell-Assisted Lipotransfer procedure.

In 2001, Dr. Patricia Zuk first published a thesis explaining how stem cells from adipose tissue were obtained and cultured successfully by their research team.

[Tissue Engineering, 7:212-228-2001]

In 2008, Dr. Yoshimura, a professor of Plastic Surgery Department of Tokyo University, proved the efficacy of a method called, 'CAL Method' (Cell-Assisted Lipotransfer) for increasing survival rate of adipose tissue. Using this technique, the survival rate of the implanted fats has increased up to 70%, where it was supposedly 30%. This method, then, was recognized as the best method of using SVF (Stromal Vascular Fraction) with simplicity and ease of usage.

[Aesthetics Plastic Surgery, 32: 48-55-2008]

Based on the theories from Dr. Zuk and Dr. Yoshimura, UNISTATION™ is specialized in SVF or ADSC isolation for CAL (Cell-Assisted Lipotransfer) therapy, which can easily be adopted to be used in private practices, general hospitals, and research centers.

ADSCs enhances Fat Survival Rate up to 70%, compared to 30% of previous Lipotransfer.





CLINICAL RESULT of CAL Method

■ Breast Augmentation

A 30 Year old woman underwent breast augmentation with CAL (310ml each breast). Her breasts were dramatically augmented with an increase in breast circumference differences **by 8cm after 24 months.**

Among 310ml injected in her breast, almost 220ml of transferred fats has survived within 24 months.

Injection amount: 310ml
Result After 2 Years: 220ml
Survival Rate: 70.9%

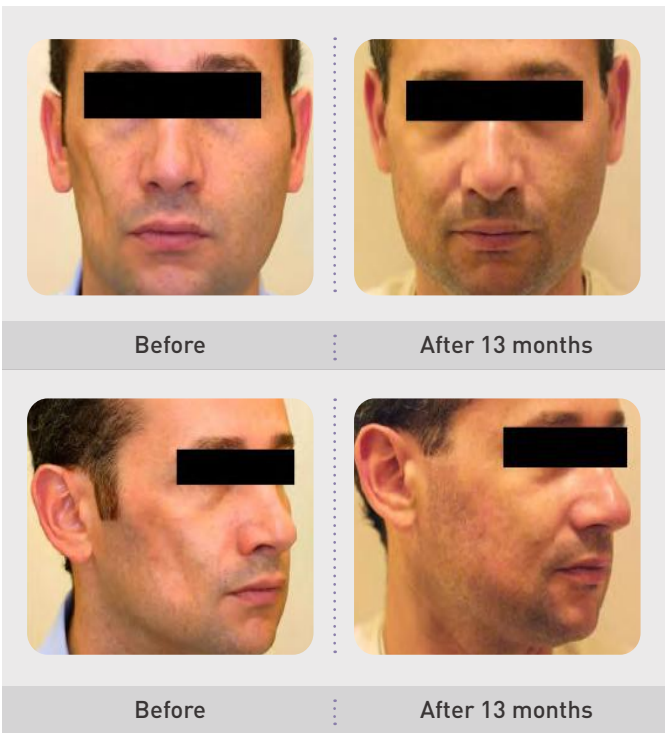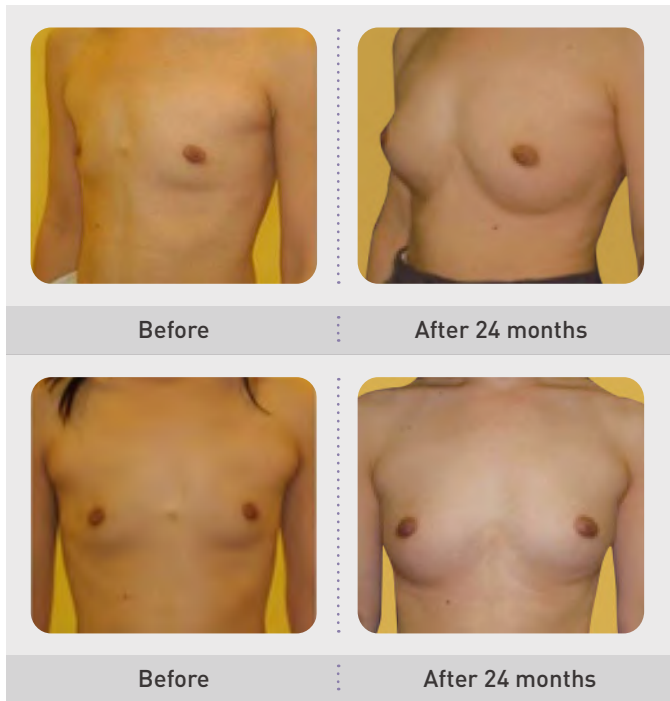

■ Lupus Profundus (a.k.a. Parry- Romberg Syndrome)

	Non CAL	Non CAL	CAL	CAL
Sex	M	F	M	F
Age	25	42	35	33
Diagnosis	LEP	LEP	PRS	LEP
Severity rating	3	5	4	5
Period	10	9	13	9
Clinical result	Fair	Good	Excellent	Good
Complications	X	Adipose Necrosis	X	X

This comparison chart clearly shows that CAL therapy resulted in a much higher outcome for clinical results. Whereas non-CAL therapy resulted in serious Adipose Necrosis.


Composition & Consumables




UNISTATION



Centrifuge

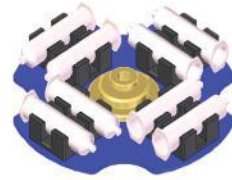

Shaker

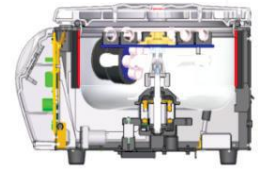

Incubator



UNIKIT (Uni-Syringe)



Transfer



Cap



Closer

ADSC Isolation Procedure




1. Preparation



2. Fat Washing



3. Collagenase Digestion

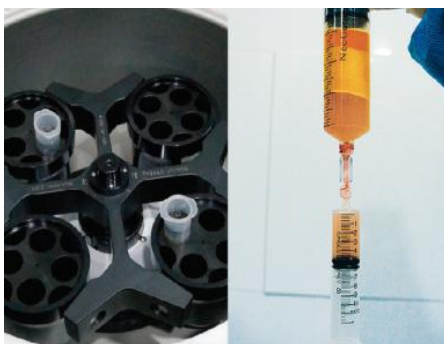

4. SVF Collection



5. Neutralization and Washing of SVF




6. Final Result



NeoGenesis PRP Series

Specially designed to extract highly concentrated PRP from blood

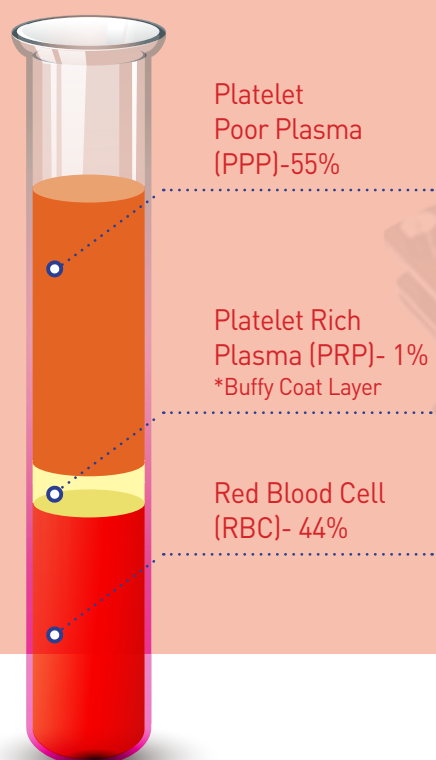

Why NeoGenesis PRP Series?

Platelets, one of the blood elements that flow along the blood stream in an inactivated form gets activated by damages of tissues and perform the original function of wound healing. As PRP therapy maximizes wound healing ability by concentrating the platelets, higher concentration of PRP are derived in a more effective PRP therapy.

Normally, general PRP therapies use about 4~6 times of concentrated PRP, but NeoGenesis PRP Series allows a higher extraction rate from the most highly concentrated PRP (maximum 48times).



Neo PRP 10cc



Genesis PRP 15cc



Pure PRP 20cc



Cure PRP 30cc



Omni PRP 60cc

Special Features of NeoGenesis PRP Series

High concentration of PRP & Easy extraction process

- The maximum concentration of PRP regardless of the operator's skill

Curved neck design to reduce the possible cell losses

- Reduction in the cell loss – The maximum ratio of PRP

PRP and PPP extraction from a single kit by one-time centrifugation

- Controls the volume of plasma according to the treatment plan

Closed system to prevent from air contamination possibility


- No need of a needle when collecting PRP or PPP

Competitive price

- The highest quality of PRP tubes along with its competitive price



Application of NeoGenesis PRP Series

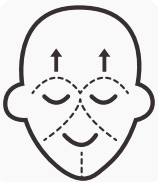

Plastic Surgery / Dermatology

- Skin rejuvenation
- Autologous fat grafting
- Acne scar regeneration



Orthopedics / Rehabilitation Medicine

- Pain relief after surgery
- Sports rehabilitation
- Joint inflammation healing
- Intractable tendon damage healing



Dentistry

- Tooth extraction & Implant
- Bone grafting



Hair Loss Prevention

- Hair loss prevention
- Hair root reinforcement



Gynecology

- Boost the endometrium cells during embryo transfer



Veterinary

- Maximizes therapeutic action after surgery
- Tendon regeneration



Neo PRP 10cc



High concentration rate of Platelets (Maximum 7~8 times)

Certificate Registration	Neo PRP Tube Size	Centrifuge Compatibility	Recommended Usage
CE, GMP, PCT ISO 13485, KFDA	Diameter: 1.6cm Height: 12.9cm Weight: 10g	Angle type Swing type	Dermatology Plastic Surgery Gynecology



Genesis PRP 15cc



High concentration rate of Platelets (Maximum 12~13 times)

Certificate Registration	Genesis PRP Tube Size	Centrifuge Compatibility	Recommended Usage
CE, GMP, PCT ISO 13485, KFDA	Diameter: 3.3cm Height: 11.1cm Weight: 24g	Swing type	Dermatology Plastic Surgery Gynecology Dentistry



Pure PRP 20cc



High concentration rate of Platelets (Maximum 15~16 times)

Certificate Registration	Pure PRP Tube Size	Centrifuge Compatibility	Recommended Usage
CE, GMP, PCT ISO 13485, KFDA	Diameter: 3.8cm Height: 12.2cm Weight: 43g	Swing type	Orthopedics Tendon Damage Rehabilitation Medicine



Cure PRP 30cc



High concentration rate of Platelets (Maximum 24~25 times)

Certificate Registration	Cure PRP Tube Size	Centrifuge Compatibility	Recommended Usage
CE, GMP, PCT ISO 13485, KFDA	Diameter: 3.8cm Height: 12.2cm Weight: 43g	Swing type	Orthopedics Osteoarthritis Intractable Tendon Damage Rehabilitation Medicine



Omni PRP 60cc



High concentration rate of Platelets (Maximum 7~8 times)

Certificate Registration	Omni PRP Tube Size	Centrifuge Compatibility	Recommended Usage
CE, GMP, PCT ISO 13485, KFDA	Diameter: 3.8cm Height: 12.7cm Weight: 45g	Swing type	Hair Implant Hair Loss Treatment

**** Omni PRP provides maximum of 45~48 times of concentration ratio, but as Omni PRP extracts 6cc of concentrated Platelets, its actual concentration ratio is about 7~8 times.**

NeoGenesis PRP Series Comparison Chart

	Neo PRP	Genesis PRP	Pure PRP	Cure PRP	Omni PRP
Required Volume of Blood	10 mL	15 mL	20 mL	30 mL	60 mL
Concentration of PRP	Maximum 7~8 times	Maximum 12~13 times	Maximum 15~16 times	Maximum 24~25 times	Maximum 7~8 times
Extracted Volume of PRP	1 mL	1 mL	1 mL	1 mL	6 mL



PRP Activator

PRP activator to maximize the effect of PRP therapy


Activating Method

PRP Activator was developed based on the idea of how Platelets become activated through narrow blood vessels in the rapid blood stream (Natural Activation).


- Uses natural activation (physical method), No chemicals included
Can prevent side effects due to non-chemical usage
- Quantitative increase of Growth Factors
Maximizes the effect of PRP Therapy
- Relatively reduces the volume of blood required compared with inactivated blood
Increases patient's convenience




Preliminary Experiment

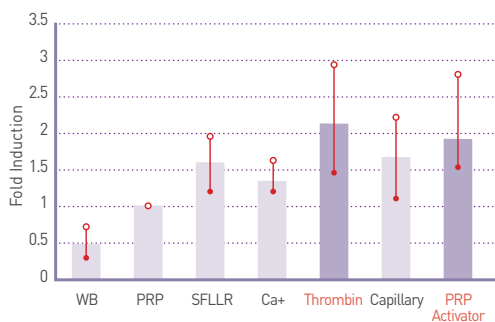

Epidermal Growth Factor

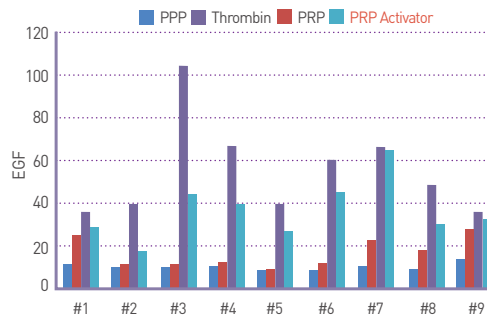

Vascular Endothelial Growth Factor

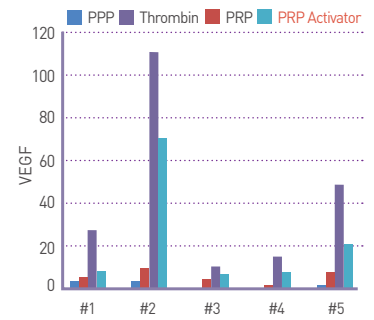

The graphs above show that the effects of PRP Activator is equivalent to Thrombin, depending on the patient's condition.

PRP Activator Instruction Manual



1 Connect two luer-lock syringes [one empty and one with PRP] to each end of PRP Activator.



2 Push the plunger back and forth.



3 Repeat Step 2 for about 20 to 30 times.



4 PRP is activated!



Genesis Centrifuge

Medical centrifuge for general laboratory applications

Advantages of Genesis Centrifuge

- **Safety**
 - CE as Class II Medical centrifuge
 - Automatic recognition and alarms for imbalance over-speed and over-heating
- **Soft Spin Key**
 - Gentle acceleration and deceleration
- **Low noise Level**
 - 52dB with unique airflow damping and elimination of sample heating
 - Double layered door construction (ABS / Steel), Decreased heat and noise transmission



Compatibility and Specification of Genesis Centrifuge

Type	Standard Edition 15cc (GRS-S-50-4)	Special Edition 30cc (HRS-60-4)
Compatible Kits	<ul style="list-style-type: none"> - Neo PRP (10cc) - Genesis PRP (15cc) - Pure PRP (20cc) - Conical Tube 15cc, 50cc 	<ul style="list-style-type: none"> - Neo PRP (10cc) - Genesis PRP (15cc) - Pure PRP (20cc) - Cure PRP (30cc) - Omni PRP (60cc) - Conical Tube 15cc, 50cc
Components	<ul style="list-style-type: none"> - Rotor: GRS-S-50-4 - Bucket: GLB-50-50 (Genesis PRP) - Adapter: GAS-U15 (50) (NEO PRP) 	<ul style="list-style-type: none"> - Rotor: HRS-6-4 - Bucket: HLB-60 (Cure PRP) - Adapter: GAS-30 (HLB 60) (GENESIS PRP) - Adapter: GAS-10 (HLB 60) (NEO PRP)

DIVISION	SPECIFICATION
MAX. RPM / RCF	4,000 RPM / 2,700 xg
MAX. Capacity	4 x 100 mL
Time Control	Timed < 100min or continuous
RPM / RCF Conversion	Yes
Noise level (dB)	± 52
ACC / DEC	± 20 / ± 25
SOFT Start / Stop	Yes
Program Memory	10
Imbalance Cutout	Yes
Safety Lid Lock	Yes
Lid Drop Protection	Yes
Automatic Door Release at Completion	Yes
Power Supply (V / Hz)	220 / 50~60 (110V optional)
Power Requirement (VA)	300
Dimension (W x D x H, mm)	375 x 480 x 260
Weight without Rotor (Kg)	19.5



Plasma Gel

Autologous filler maker from Platelet Poor Plasma (PPP, or Plasma) for contouring effect.

Advantages of Plasma Gel

- **Safety**

PPP is extracted from the patient's own blood

- **Combination therapy with PRP**

Can be easily combined with PRP therapy to strengthen the efficacy of procedure

- **Quick & Simple Procedure**

One button protocol (automatic process)

Total preparation time: approximately 18 minutes

- 70 degrees Celsius → 5 minutes

- 90 degrees Celsius → 8 minutes

- Cooling (approximately 5 minutes) / Natural Cooling



Combined Therapy and Specification of Plasma Gel



When activated PRP is added into the Plasma Gel, duration of the self-made filler increases, while the growth factors contributing to the rejuvenation of the skin tissues improve.

DIVISION	SPECIFICATION
Dimension	18(L) x 16(W) x 15(H) cm
Weight	1.9 kg
Production Capacity	6cc at once
Syringe Holder Capacity	1cc syringe x 6 holes
Display Indicating Current Temperature	Yes
Power/Emergency Stop Buttons	Yes



REGENOVUE

HA filler series to correct volume loss and help create a more youthful appearance

What is REGENOVUE?

REGENOVUE series consists of mono-phasic Hyaluronic acid fillers (the safest material existing in the human body) that is 100% cross-linked for a longer lasting duration. Regenovue is a all natural, safe and effective filler that offers a full line up of 6 models (with/ without Lidocaine), for diverse treatment areas depending on the molecular size and lidocaine contents.



REGENOVUE FINE



REGENOVUE FINE PLUS




REGENOVUE DEEP



REGENOVUE DEEP PLUS



REGENOVUE SUB-Q

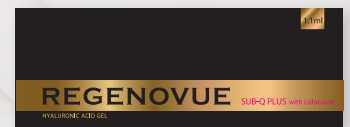

REGENOVUE SUB-Q PLUS

Features of REGENOVUE

High Quality Hyaluronic Acid

- Non-animal base
- Bacterial extract natural substance
- High Purity (Endotoxin <0.0015IU/mg>
- Biodegradable

Highly Cross-linked

- Linking material:BDDE
- Specialized manufacturing technique
- Increases volume and duration

High Viscosity

- Long-lasting duration
- Supports skin texture stably to maintain elasticity

Full Line-up Models

- Various choice and accurate application
- Includes Plus series for lidocaine





Comparison of Structure


REGENOVUE

The even molecule size of Structure of REGENOVUE indicates a specialized manufacturing technique of NeoGenesis Co., Ltd.

This will increase the duration and give smooth injecting feeling to the doctors while reducing patient's discomfort.

Structure of J Product

The picture indicates that the J product's molecular sizes are comparatively less equal and less crosslinked than REGENOVUE.

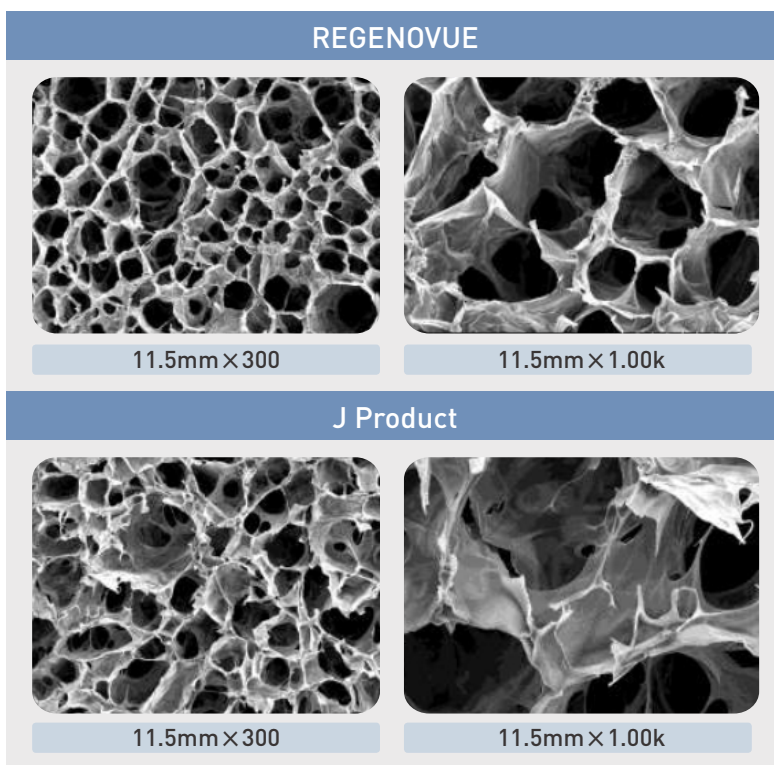

REGENOVUE Product Detail

Plus models contain lidocaine

	REGENOVUE Fine / Fine Plus	REGENOVUE Deep / Deep Plus	REGENOVUE Sub-Q / Sub-Q Plus
Ingredient	HA 24mg/ml	HA 24mg/ml	HA 24mg/ml
Application	Superficial wrinkles, periorbital lines, glabellar lines, neck folds	Middle/deep wrinkles, nasolabial folds, face, lips, chin, forehead	Deep wrinkles, nose surgery
Injecting Area	Upper Dermis	Lower Dermis, Upper Sub-Q	Lower Sub-Q
Syringe	1.0ml / 1.1ml	1.0ml / 1.1ml	1.0ml / 1.1ml
Needle Size	30G	27G	27G
Duration	8~12 months	12~18 months	12~18 months
Storage Temp.	2~25°C	2~25°C	2~25°C
Shelf Life	24 months	24 months	24 months
Packing	1 syringe / 1 x 30G, 1 x 27G needles	1 syringe / 2 x 27G needles	1 syringe / 1 x 27G, 1 x 25G needles
Lidocaine	X / O	X / O	X / O



REGENOVUE Aqua Shine

HA Mesotherapy for a more concentrated and rejuvenated hydrated appearance

What is REGENOVUE Aqua Shine?

REGENOVUE Aqua Shine Treatment is an anti-aging treatment to rejuvenate the skin by injecting HA evenly and directly into the dermis, resulting in brighter and moister skin by concentrating moisture in skin as HA has the effect of attracting moisture.



REGENOVUE Aqua Shine



REGENOVUE Aqua Shine Plus

Effect of REGENOVUE Aqua Shine

Fine Wrinkle Removal

Increase Skin Elasticity

Moisturizing Effect

Anti-Aging Care

Improve Skin Tone

Immediate Result

Pore Tightening

REGENOVUE Aqua Shine Specification

Name	REGENOVUE Aqua Shine
Features	Pure HA
Advantage	Immediate Effect
Main Composition	Hyaluronic Acid 15mg/ml
Specification	3.0ml × 3ea

Name	REGENOVUE Aqua Shine Plus
Features	Cross-linked HA
Advantage	Long-lasting Effect
Main Composition	Cross-linked Hyaluronic Acid 15mg/ml
Specification	3.0ml × 3ea



Neo Filler Cannula

Specially designed convenient filler cannulas for the most effective results

Advantages of Neo Filler Cannula

- Smooth Injection
- Smooth Surface of Cannula
- Wider Inner Diameter of Cannula
- Diverse Types (14g – 30g & 25mm – 200mm)
- KFDA Certified
- GMP Certified
- ISO 13485 Certified
- CE Applied



Neo Filler Cannula Model List

Name	Cannula Gauge	Cannula Length	Model Name	Remarks
NEO Filler Cannula	30G	25mm	MC-30025	1. MOQ - 50 pcs (1 Box) for each model 2. Marked Cannula can be applicable if specially requested
		30mm	MC-30030	
	27G	30mm	MC-27030	
		40mm	MC-27040	
	25G	40mm	MC-25040	
		50mm	MC-25050	
	24G	50mm	MC-24050	
		30mm	MC-23030	
	23G	50mm	MC-23050	
		70mm	MC-23070	
NEO Fat Graft Cannula	22G	50mm	MC-22050	
		50mm	MC-21050	
	20G	50mm	MC-20050	
		70mm	MC-20070	
	19G	50mm	MC-19050	
		70mm	MC-19070	
	18G	50mm	MC-18050	
		70mm	MC-18070	
NEO Breast Cannula	17G	70mm	MC-17070	1. MOQ - 20 pcs (1 Box) for each model 2. Marked cannulas can be applicable if specially requested
		100mm	MC-17100	
		120mm	MC-17120	
	14G	150mm	MC-17150	
		200mm	MC-14200	1. MOQ - 10 pcs (1 Box) for each model



Neo Thread Series

Specially designed lifting PDO (Polydioxanone) threads with wide variety of models for outstanding results

What is Neo Thread Series?


Neo Thread Series are widely diverse as it can be used not only for the face, but also for the body such as breast, abdomen, legs, buttocks, and more as well as treatment.



Advantages of Neo Thread Series

Neo Thread Series includes a wide variety of threads to provide the most satisfying result and beauty. **Neo Thread Series** is made from Polydioxanone material (PDO) with types of series that includes:

- Neo Mono Thread – Mono and Screw type threads
- Neo Multi Thread – Combination with Cell therapies
- Neo Cog Thread – Available in multi-directional barbed thread, and cannula
- Neo Combi Thread – Available in multi-directional barbed thread, and cannula
- Neo Nose Thread – Specially designed for nose with cannulas
- Neo Eye Thread Cannula – Specially designed for eye with cannulas
- Double Needle Thread – For professional modifications
- Molding Thread Cannula – Specially designed thread cannula pressed with heat to mold barbed cog structure for face and body

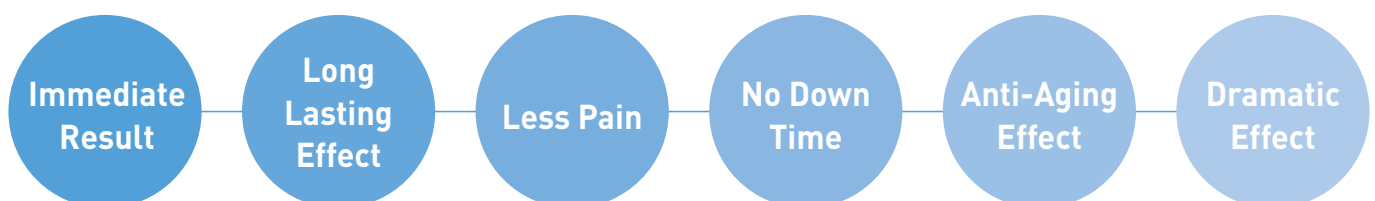






Application & Combination Therapy

Face Lifting by Combi + Multi Thread

Combi for entire face contour, nasolabial line and double chin

Multi for deep and fine wrinkles on forehead, for thin and sensitive skin around eyes and lips

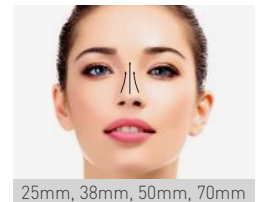
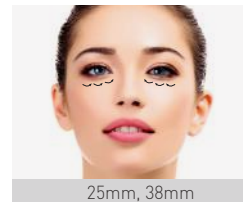
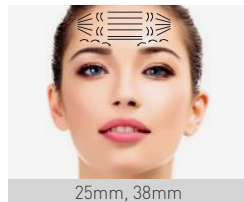


Face Lifting by Cog + Combi + Multi Thread

Cog for entire face contour and nasolabial line

Combi for nasolabial line and double chin

Multi for deep and fine wrinkles

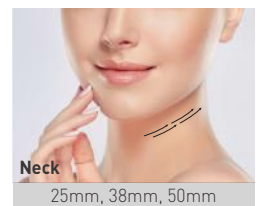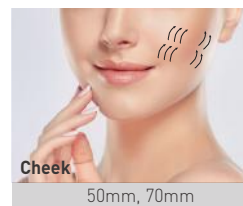
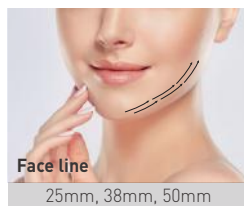


Body Contour Lifting by Combi Thread

Combi for breast lifting

Combi for abdominal lifting

Combi for buttock lifting

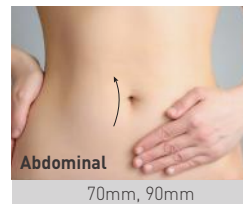
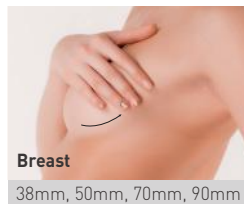


Specially designed Nose, Eye Cannulas, and Double Needle Threads

Cannula for nasal ridge, and nasal base

Cannula around the eyes



Multi-directional coggled thread for total lifting corrections

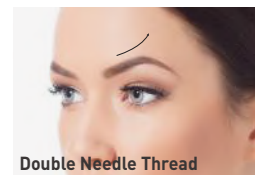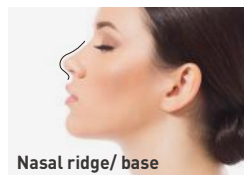

Neo Mono Thread

For face, double chin, body contouring (neck, breast, buttock), etc

NEO MONO THREAD is made up of PDO material thread that induces a self-healing process by inserting thread into the skin. Neo Mono Thread comes in two types of thread, Mono and Mono Screw, for collagen production, anti-wrinkle, and skin rejuvenation effects.

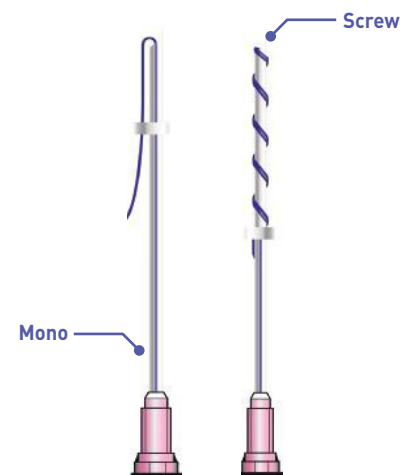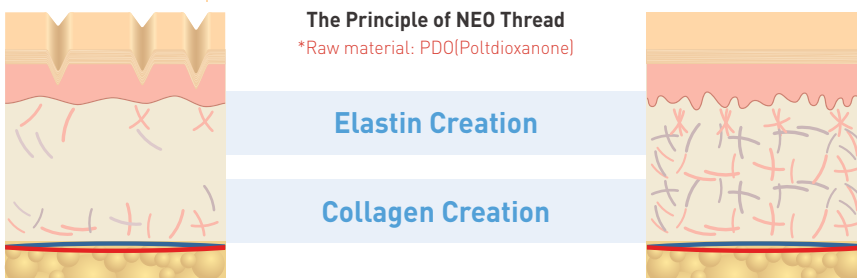

Neo Mono and Screw Thread

Neo Multi Threads

For stronger lifting with scaffolding effect in diverse cell therapies for face and body

Multiple threads combined to a single thread to stimulate collagen and cell regeneration and to fill up hollow cheek, chin, eyes, as well as removing deep and fine wrinkles. NEO MULTI THREAD can produce more collagen and lifting to the skin with combined cell therapy due to its engineered scaffolding structure.

Cell Therapy

- Combined cell therapy with Multi Thread by injecting stem cells and/or PRP (Platelet Rich Plasma)
- Higher probability (possibility) of cell viability and engraftment (Multiple threads act as Scaffolds)



Collagen with Mono Thread



Collagen with Multi Thread



NEO Multi Thread

Neo Cog Thread


For an all-time solution to uplift face and body with a thorn-shaped cog thread

NEO COG THREAD provides the best way to lift sagging cheek, chin, jaw line and body contouring through its most simple procedure. Neo Cog Thread, by offering two types of directions of Cog thread will give you excellent results in making deep wrinkles and furrows uplifted.

- Raw material : PDO(Polydioxanone)
- Diameter : USP 0-0 / 1 / 2
- Shape : Thorn-shaped Cog
- Direction : Spiral / Bi-Directional



Model Name	Needle	Thread Length	Thread USP
C19100	19G 100mm	140mm	0-0
C19100-C (U-Type Cannula)	19G 100mm	140mm	0-0
C19100B-C (Bi-Directional/U-Type Cannula)	19G 100mm	140mm	0
C18100-C (U-Type Cannula)	18G 100mm	150mm	2

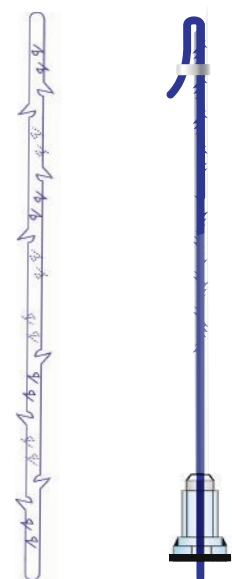

NEO Cog Thread
* An illustrated image of 19G




Neo Combi Thread

For face, double chin, body contour and lifting(neck, breast, buttock), etc


NEO COMBI THREAD includes deep cogs with thickness that provides stronger pulling force for effective result in lifting sagging skin. Minimal quantity of Combi threads are required to make amazing result through its multi-directional cog thread.

- Raw material : PDO(Polydioxanone)
- Diameter : USP 2-0 / USP 3-0
- Shape : Thorn-shaped Cog
- Direction : Multi-Directional

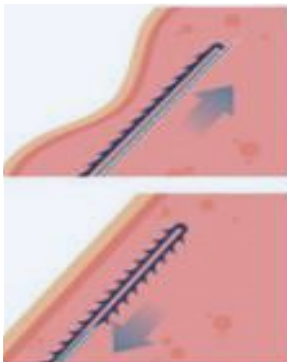

The Principle of NEO Combi Thread

Model Name	Needle	Thread Length	Thread USP
B2360	23G 60mm	90mm	3-0
B2390	23G 90mm	150mm	3-0
B2160	21G 60mm	90mm	2-0
B2190	21G 90mm	150mm	2-0
B2360-C (U-Type Cannula)	23G 60mm	90mm	3-0
B2390-C (U-Type Cannula)	23G 90mm	150mm	3-0
B2160-C (U-Type Cannula)	21G 60mm	90mm	2-0
B2190-C (U-Type Cannula)	21G 90mm	150mm	2-0

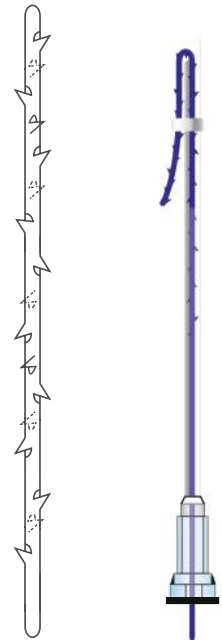

NEO Combi Thread


* An illustrated image of 21G/23G

Neo Nose Thread


Specially designed L-type cannula for nasal ridge, and nasal base lifting through an anchoring effect by the thorn-shaped cog thread

NEO NOSE THREAD is a thorn-shaped cog thread that stimulates a lifting effect by anchoring to the dermis of the nose when inserted. Neo Nose Thread offers lifting of the nasal ridge, and shrinkage of the nasal base by re-shaping the nose. A safe and simple treatment compared to Rhinoplasty.

* L-Type Blunt Cannula needles are used for nose thread models

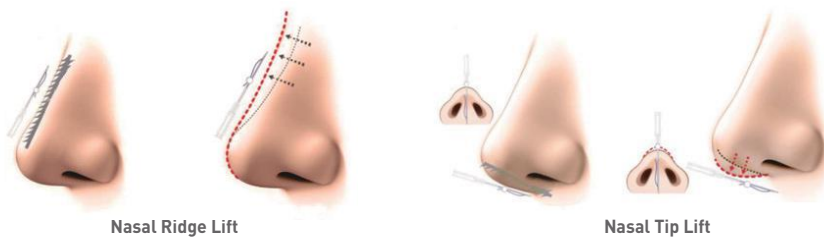

Model Name	Needle	Thread Length	Thread USP
C1938-C (L-Type Cannula)	19G 38mm	50mm	2
C1950-C (L-Type Cannula)	19G 50mm	80mm	2

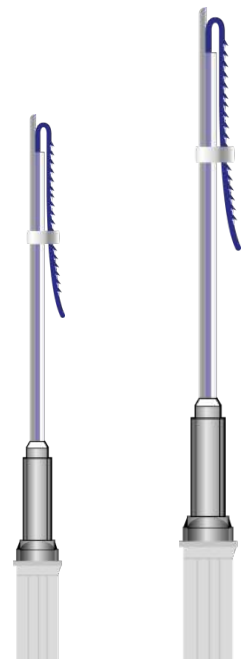

NEO Nose Thread

* An illustrated image of C1950-C

Neo Eye Thread Cannula

Specially designed U-type Eye Thread Cannulas to remove wrinkles and bags around the eyes

NEO EYE THREAD CANNULA reduces the sagging around the eyes while minimizing bleeding, swelling, and bruising. The Mono and Screw type Eye cannulas continuously gives increased production of collagen.



Model Name	Needle	Thread Length	Thread USP
M2938-C (U-Type Cannula)	29G 38mm	50mm	6-0
S2938-C (U-Type Cannula)	29G 38mm	50mm	6-0



NEO Eye Thread Cannula


*Mono and Screw Eye cannulas in 29G

Neo Double Needle Thread

A Specially designed multi-directional cogged thread for total lifting corrections

NEO DOUBLE NEEDLE THREADS insert without the need of a cannula, with both needles inserted at the midpoint of the path and driven upwards to their respective ends. Standard techniques can be used for forehead/brow lifts, cheek /neck lifts, and jaw lines correction.

Model Name	Needle	Thread Length	Thread USP
D19110S (Short)	19G 110mm	290mm	0-0
D19110L (Long)	19G 110mm	440mm	0-0

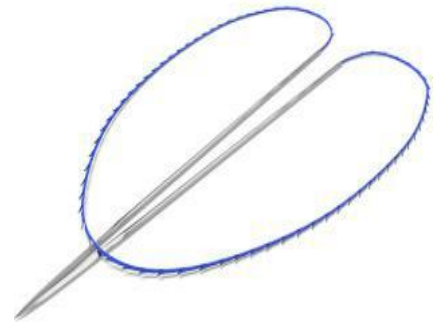

Double Needle Thread

*Illustrated Image of Double Needle Thread in 19G

Neo Molding Thread Cannula

Specially designed molding cog thread cannula for face contouring (sagging cheek & chin), and body contouring (neck, breast, and buttock)

NEO MOLDING THREAD CANNULA offers much more advantages in terms of lifting, strength, duration, and cost. As it is pressed with heat to mold its barbed cog structure, Molding Thread Cannula maintains its original form of thickness intact with barbs.

Molding Threads have been developed by a mold-pressing heat process, which conserves the integrity of the thread. Before the molding process, the threads have large thick "bumps" along the threads. As these bumps are pressed, the original large and thickness of the thread do not change, and shape these bumps into barbs.

Model Name	Needle	Thread Length	Thread USP
MC18100-C (L-Type Cannula)	18G 100mm	150mm	3



Molding Thread Cannula

*Illustrated Image of Molding Thread Cannula in 18G



Category List for Neo Thread Series

Name	Needle Size	Needle Length	Thread Length	Thread USP	Model Name	Remarks
NEO Thread (Mono)	30G	25.4mm	30mm	7-0	M3025	Mono - Mono Filament / PDO
	30G	38mm	50mm	7-0	M3038	
	29G	25.4mm	30mm	6-0	M2925	
	29G	38mm	50mm	6-0	M2938	
	29G	50mm	70mm	6-0	M2950	
	29G	60mm	90mm	6-0	M2960	
	27G	38mm	50mm	5-0	M2738	
	27G	50mm	70mm	5-0	M2750	
	27G	60mm	90mm	5-0	M2760	
	26G	38mm	50mm	5-0	M2638	
	26G	50mm	70mm	5-0	M2650	
	26G	60mm	90mm	5-0	M2660	
	26G	90mm	150mm	5-0	M2690	
NEO Thread (Screw)	30G	25.4mm	30mm	7-0	S3025	Screw - Mono Filament / PDO
	30G	38mm	50mm	7-0	S3038	
	29G	25.4mm	30mm	6-0	S2925	
	29G	38mm	50mm	6-0	S2938	
	29G	50mm	70mm	6-0	S2950	
	27G	50mm	70mm	5-0	S2750	
	27G	60mm	90mm	5-0	S2760	
	26G	60mm	90mm	5-0	S2660	
	26G	90mm	150mm	5-0	S2690	

Name	Needle Size	Needle Length	Thread Length	Thread USP	Model Name	Remarks
NEO Multi Thread (Twin)	27G	38mm	50mm	(7-0) X 2	T2738	Twin - 2 Twisted Filaments / PDO
	27G	50mm	70mm	(7-0) X 2	T2750	
	27G	60mm	90mm	(7-0) X 2	T2760	
	26G	38mm	50mm	(7-0) + (6-0)	T2638	
	26G	50mm	70mm	(7-0) + (6-0)	T2650	
	26G	60mm	90mm	(7-0) + (6-0)	T2660	
	26G	90mm	150mm	(7-0) + (6-0)	T2690	
	25G	60mm	90mm	(6-0) X 2	T2560	
	25G	90mm	150mm	(6-0) X 2	T2590	

Name	Needle Size	Needle Length	Thread Length	Thread USP	Model Name	Remarks
NEO Multi Thread (Triple)	26G	38mm	50mm	(7-0) X 3	TR2638	Triple - 3 Twisted Filaments / PDO
	26G	50mm	70mm	(7-0) X 3	TR2650	
	26G	60mm	90mm	(7-0) X 3	TR2660	
	26G	90mm	150mm	(7-0) X 3	TR2690	
	25G	60mm	90mm	(7-0) X 2 + (6-0)	TR2560	
	25G	90mm	150mm	(7-0) X 2 + (6-0)	TR2590	
	23G	60mm	90mm	(6-0) X 3	TR2360	
	23G	90mm	150mm	(6-0) X 3	TR2390	
NEO Multi Thread (Tornado Screw)	26G	38mm	50 mm	(7-0) X 2	TS2638	Tornado Screw - 2 Twisted Filaments / PDO
	26G	50mm	70 mm	(7-0) X 2	TS2650	
	26G	60mm	90 mm	(7-0) X 2	TS2660	
	26G	90mm	150 mm	(7-0) X 2	TS2690	
	25G	60mm	90 mm	(6-0) X 2	TS2560	
	25G	90mm	150 mm	(6-0) X 2	TS2590	

Name	Needle Size	Needle Length	Thread Length	Thread USP	Model Name	Remarks
Neo Eye Thread (Screw / Cannula)	29G	38mm	50mm	6-0	M2938-C	Mono - Mono Filament / U-Type Cannula / PDO
	29G	38mm	50mm	6-0	S2938-C	Screw - Mono Filament / U-Type Cannula / PDO
NEO Cog Thread (360 Degree / Cannula)	19G	70mm	110 mm	0-0	C1970	Cog - 360 Degree Spiral Cogged Filament / PDO
	19G	100mm	140 mm	0-0	C19100	
	19G	70mm	110mm	0-0	C1970-C	Cog - 360 Degree Spiral Cogged Filament / PDO / U-Type Cannula
	19G	100mm	140mm	0-0	C19100-C	
NEO Cog Thread (Bi-Directional / Cannula)	19G	100mm	140mm	1	C19100B-C	Cog - Bi-Directional Cogged Filament / PDO / U-Type Cannula
	18G	100mm	140mm	2	C18100-C	
NEO Nose Thread	19G	38mm	50mm	2	C1938-C	Barbed / PDO / L-Type Cannula
	19G	50mm	80mm	2	C1950-C	
NEO Combi Thread (Multi-Directional)	23G	60mm	90mm	3-0	B2360	Combi - Multi-Directional Cogged Filament / PDO
	23G	90mm	150mm	3-0	B2390	
	21G	60mm	90mm	2-0	B2160	
	21G	90mm	150mm	2-0	B2190	
	23G	60mm	90mm	3-0	B2360-C	Combi - Multi-Directional Cogged Filament / PDO / U-Type Cannula
	23G	90mm	150mm	3-0	B2390-C	
	21G	60mm	90mm	2-0	B2160-C	
	21G	90mm	150mm	2-0	B2190-C	
NEO Double Needle Thread (Multi-Directional)	19G	110mm	290mm	0-0	D19110S	Double Needle - Multi-Directional Cogged Filament / PDO
	19G	110mm	440mm	0-0	D19110L	
NEO Molding Thread Cannula [Bi-Directional]	18G	100mm	150mm	3	MC18100-C	Molding Thread Cannula - Cogged Filament / PDO / L-Type Cannula



Beta Series


Beta Serum / Beta Scaffold Mask

Intensive Hydration with 100% Activate Beta Glucan



■ ■ What is Beta Glucan, the main ingredient of Beta Series?

- Natural polymer produced by natural botanical resources
- Well-known to boost the immune system, reconstruct inner and outer surface of the skin, prevent skin infections, and irritation (anti-inflammatory effects)
- Moisturizes the skin 20% more than Hyaluronic Acid (HA)
- Beta Glucan rapidly reaches to the dermis layer without any injections

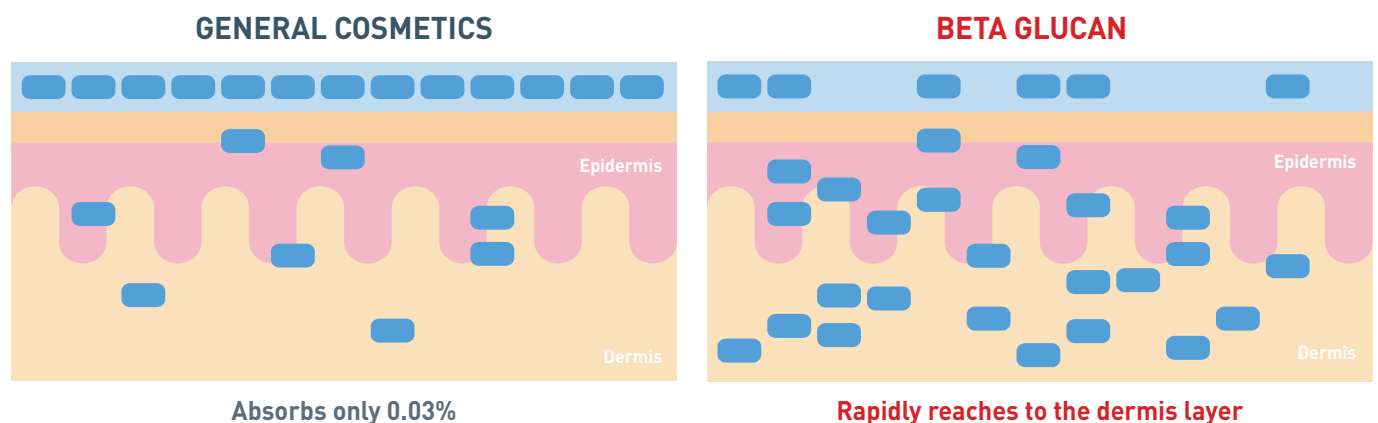

[Meyer, Wilfried MD, Schonnagel, Britta MD. 2006, A note on integumental [1→3][1→6][Beta]-D-glucan permeation, using the porcine ear skin model. Journal of Cosmetic Dermatology. 5(2):130-134, 2006]

Beta Scaffold Mask

- Provides a great cooling effect which will be very helpful for calming burns and lacerations
- Prevents swelling and redness
- Recommended post-op care for various dermatological operations
- For dry and inflamed skin for home care



Facial gel mask sheets × 10 ea per box

Beta Serum

- Moisturizes the skin by reaching to the dermis layer without injection
- Provides 20% more moisturizing effect than HA
- Tightens the pore and controls overall skin tone
- Recommended post-op care for various dermatological operations
- For dry and inflamed skin for home care on a daily basis



Beta Serum Bottles (15ml) × 3

(Male, 23)




Before (2010. 08. 07)



After 3 months (2010. 11. 20)

■ The Effect of Beta Series on Severe Scalp Psoriasis

Beta Serum provides skin enhancement especially for severe psoriasis due to anti-inflammatory and moisturizing effects, when used regularly. The pictures on the side clearly show that the patient suffered from severe psoriasis. By applying Beta Serum for 3 months, the patient had a satisfying result.

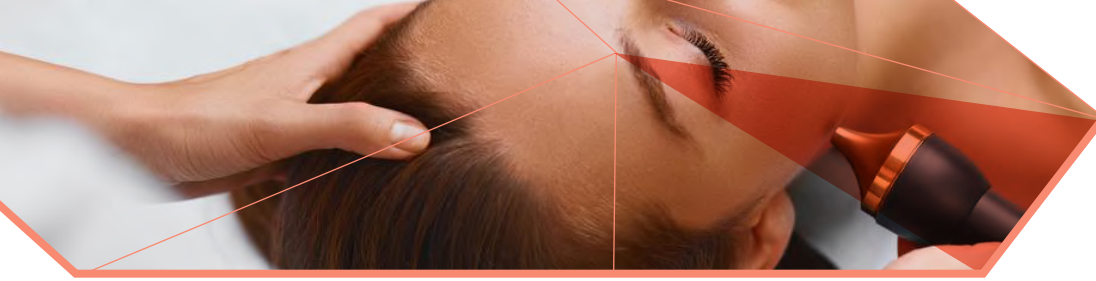

SPOTLESS



Radio Frequency device to erase blemishes

What is Spotless?

- Removes small blemishes on your skin with RF (Radio Frequency)
- Minimal-invasive, painless, no downtime
- Immediate and permanent effect with one-time treatment
- Minimal damage on skin surface
- Reasonable cost compared to other expensive and heavy laser equipments
- CE, ISO 13485, GMP Certified
- Non-chargeable Device




Application Area

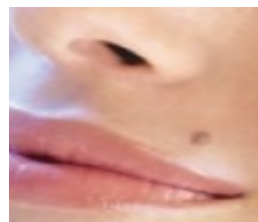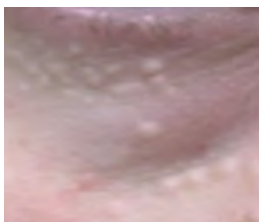

- Moles
- Freckles
- Warts
- Milium
- Syringoma

Directions

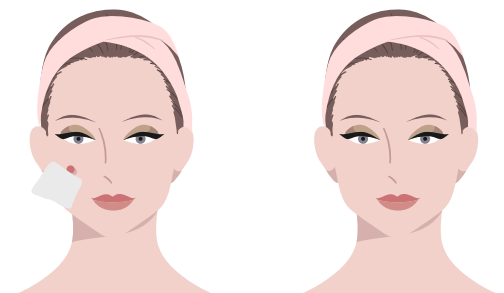

1. Wipe treating area with alcohol gauze.
2. Apply Neo-Cain (lidocaine) cream to treating area.
3. Control the power and frequency to best fit your application area.
4. Apply the device with the right balance and aim on to the skin.

Product Details




Package Box



Main Device



Adaptor

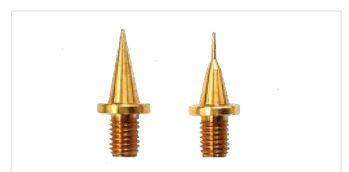

Tips

RIOX

CO₂ Carboxyl facial mask to add oxygen and nutrients to your face

What is RIOX?

RIOX is a facial mask, made to utilize Carboxyl therapy technique, a professional medical technique that has been performed only by dermatologists, plastic surgeons, and aesthetic clinics.

Its unique skin therapy technique through CO₂ provides the skin with rich nutrients of oxygen, detoxification, and balanced skin tone. In addition, restores hydration and adds brightness to the skin.





Effects of RIOX

- Improves and restores skin brightness
- Detoxification of dermis layer
- Balances oil & moisture level
- Pore tightening
- Hydration and moisturizing of the skin
- Balanced skin tone
- Reduction of pigmentation, freckles, spots
- Clearing of acnes and post-acne marks

What is Bohr Effect?

Bohr Effect is increase in CO₂ that results to a decrease in pH level. Through this effect, it releases large amount of oxygen to our skin. The CO₂ will excrete waste matters from the skin and provide nutrients and oxygen. The high density of CO₂ stimulates the metabolism of the skin so that the skin will restore moisture level, skin brightness, and healthiness that is long-lasting.

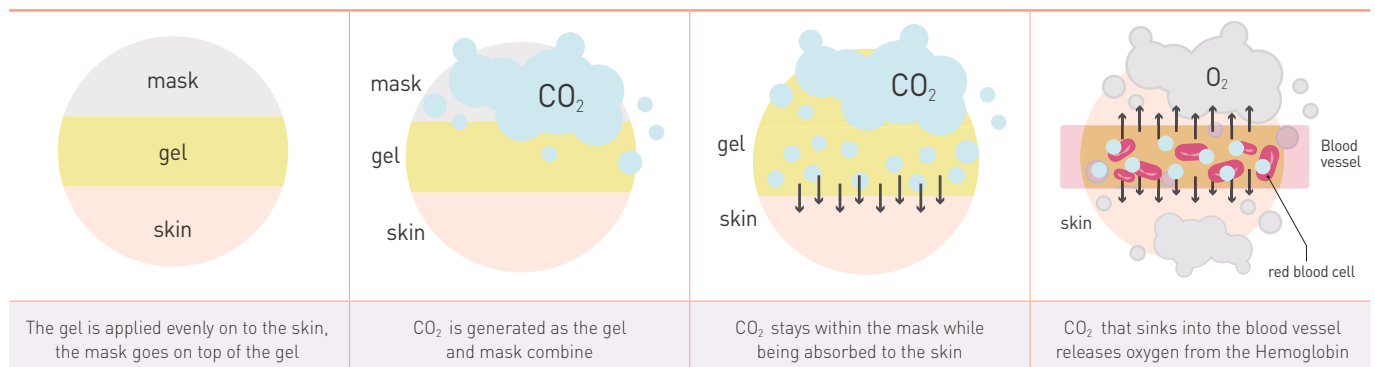

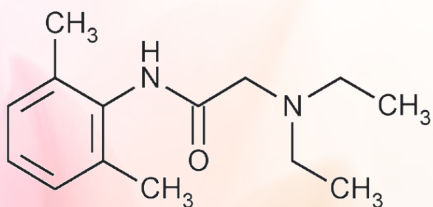

Neo-Cain Cream

Lidocaine Cream 10.56%



What is Neo-Cain Cream?

- Numbing anesthetic cream
- For loss of feeling to the skin and surrounding tissues
- To prevent and to numb pain from aesthetic procedures



Neo-Cain Cream prevents the skin to feel pain during petit procedures, while allowing patients to feel comfortable during treatment.

The following application areas include: tattooing, incisions, various skin treatments (laser, lifting thread procedures), beauty petit procedures (filler injections, moles removals), and more.

Application Area

DERMATOLOGY	UROLOGY	PLASTIC SURGERY	TATTOO

Porduct Details

Composition

1g contains

Active ingredient : Lidocaine ···· 106.5mg

Additives (Preservatives) : Methyl Paraben (KP) ···· 1.5mg

Appearance

White cream

Indications

To reduce skin sensitivity

Packing units

1 Neo Cain Tube Type (30g)

1 Neo Cain Jar Type (500g)





Address #302, 303 E&C Dream Tower, 146, Seonyu-ro, Yeongdeungpo-gu, Seoul, South Korea(07255)

Tel +82-2-583-1348 **Fax** +82-2-2038-0848 **E-mail** info@neogenesis.co.kr

Homepage www.neogenesis.co.kr **Facebook** www.facebook.com/neogenesis.coltd/